



VINNERLOVE RINSE-FREE SHAMPOO CAP (VLSHMC)

VINNERLOVE 免冲洗洗发浴帽

DIGNIFIED HAIR CLEANSING FOR CARE RECIPIENTS.
EASIER BEDSIDE CARE FOR CAREGIVERS.

让被护理者清洁更体面，让床旁护理更轻松。



CATALOGUE

目录

- PRODUCT SPECIFICATIONS	3
- KEY FEATURES	4
- COMPARED WITH TRADITIONAL CARE	5
- USAGE SCENARIOS & TARGET USERS	7
- DIRECTIONS FOR USE	8
- OTHER VINNER PRODUCT	9
- 产品参数	3
- 关键亮点	4
- 与传统床旁护理对比	5
- 使用场景与对象	7
- 使用说明	8
- 其他维尼产品	9

PRODUCT SPECIFICATIONS

产品参数

EN		CN	
Item	Specification	项目	规格
Product Name	VINNERLOVE Rinse-Free Shampoo Cap	产品名称	VINNERLOVE免冲洗洗发浴帽
Product Code	VLSHMC1	产品型号	VLSHMC1
Cap Diameter	30 CM	浴帽直径	30 CM
Pack Size	1 Cap / Pack	包装规格	1浴帽/包
Fabric Weight	Ultra-Thick Multi-Layer Nonwoven	布克重	超厚多层无纺布
Formula Type	No-rinse cleansing	配方类型	免冲洗清洁护理液
Heating Method	Microwaveable	加热方式	可微波加热
Application	Rinse-free hair cleansing	使用用途	免冲洗清洁护理头发
Skin Compatibility	Dermatologically tested	皮肤相容性	经皮肤科医生测试

KEY FEATURES

关键亮点



85GSM Ultra-Thick Multi-Layer Nonwoven Fabric

厚实多层，柔韧耐用

Made with 85GSM ultra-thick multi-layer nonwoven fabric, offering a soft yet durable structure for comfortable hair and scalp cleansing.

采用 85GSM 超厚多层无纺布，柔软舒适且厚实耐用，适合头发与头皮清洁护理。



Richly Moistened for Effective Cleansing

含液充足，清洁充分

Richly saturated with cleansing solution to help thoroughly moisten the hair and scalp, supporting a more complete and comfortable cleansing experience.

充分浸润清洁液，有助于湿润头发与头皮，使清洁过程更加充分、舒适。



Dignified Hair Care Experience

洁净体面，护理尊严

Helps care recipients maintain clean, fresh hair when conventional hair washing is difficult, supporting personal hygiene, comfort, and dignity.

在传统洗头不便时，帮助被护理者保持头发清洁清爽，提升个人卫生、舒适感与体面感。



Reduced Caregiver Burden

简化流程，减轻负担

Simplifies bedside hair cleansing by reducing the need for water, basins, towels, patient repositioning, and post-care clean-up.

简化床旁洗发流程，减少取水、脸盆、毛巾、体位调整及后续清理等步骤，减轻护理人员负担。



Soft and Comfortable Cap Fit

柔软贴合，舒适包裹

The soft cap design gently fits around the head, helping the cleansing solution reach the hair and scalp evenly while improving user comfort.

柔软浴帽设计可自然贴合头部，帮助清洁液均匀接触头发与头皮，同时提升使用舒适度。



Microwaveable for Warm Care

微波加热，温暖护理

Can be warmed in a microwave before use to provide a soothing, warm cleansing experience, especially suitable for bedside and assisted care.

使用前可通过微波炉适当加热，带来温暖舒缓的清洁体验，尤其适合床旁护理和辅助护理场景。



Dermatologically Tested

经皮肤科医生测试

Dermatologically tested with a skin-friendly formula

经皮肤科医生测试，温和亲肤。

COMPARED WITH TRADITIONAL CARE

与传统床旁护理对比

Traditional Hair Washing Method 传统护理方案

01 Prepare shampoo, towels, water basin, waterproof pads, and other supplies
准备洗发水、毛巾、水盆、防水垫等工具

02 Reposition the care recipient and maintain a fixed posture
调整受护理者体位，保持固定姿势

03 Wet the hair → Apply shampoo → Massage and cleanse → Rinse → Dry
打湿→涂抹洗发水→揉洗→冲洗→擦干

04 Check for spilled water and clean the surrounding area
检查溅出的水，清理现场

VS


VINNERLOVE Shampoo Cap VINNERLOVE免洗浴帽

01 Microwave heat if a warm and comfortable experience is desired
微波加热（若希望温热的体验）

02 Put on the shampoo cap, gently massage the scalp, then remove the cap
取出浴帽，按摩头部，取下

03 Simply towel dry if needed
用干毛巾简单擦拭

04 Minimal to no cleanup required
基本无需清理现场

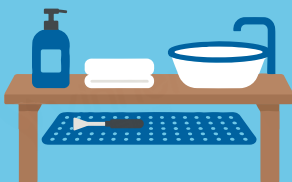

V





V



V



V

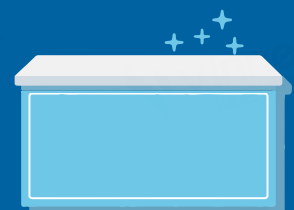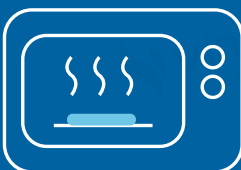

CN	方案	准备工作	洗发过程	清理工作	受护理者需配合	受护理者体验	护理风险
	VINNERLOVE	开包即用： 解放护理人员	戴上浴帽并揉搓	使用后直接丢弃	可在卧姿或坐姿下完成洗发	有安全感、 尊严感	减少飞溅风险、 减少用水、 减少共用器具 的风险
	传统护理	准备众多工具： 工作繁琐	打湿 ↓ 涂洗发水 ↓ 揉洗 ↓ 冲洗 ↓ 擦干	倒水 ↓ 清洁水盆 ↓ 处理湿毛巾和床铺	需调整头部、 颈部、 俯身或抬头	易紧张、 怕流进眼耳、 弄湿身体或床	水花飞溅、 滑倒、 污染水盆、 交叉使用 导致院内感染等
EN	Solution	Preparation	Process	Clean-up	Care Recipient Involvement	Comfort & Dignity	Risk
	VINNERLOVE	Ready to use straight from the pack; reduces caregiver workload	Simply place the cap on the head and gently massage	Dispose of after use	Hair washing can be completed while lying down or sitting	Greater comfort, dignity, and sense of security	Reduces splashing, lowers water usage, and minimizes the risk of shared equipment contamination
	Traditional Care	Requires multiple tools and complicated preparation	Wet hair ↓ Apply shampoo ↓ Scrub ↓ Rinse ↓ Dry	Empty water ↓ Clean basins, ↓ Handle wet towels and bed linens	Requires repositioning of the head and neck, bending forward or lifting the head	Can cause discomfort, anxiety, water entering eyes/ears, and wet clothing or bedding	Risks of splashing, slipping, contaminated basins, and cross-infection in healthcare settings

USAGE SCENARIOS

使用场景

EN		CN	
Scenario	Target Users	场景	适用对象
Bedside Bathing	Bedridden, post-operative, ICU and ward patients	床旁沐浴	卧床、术后、ICU及普通病房患者
Elderly Care	Seniors and people with limited mobility	老人护理	老年人及行动不便者
Long-Term Care	Nursing home and care facility residents	长期照护	护理院及照护机构住民
Home Care	Family caregivers and home recovery patients	家庭护理	家庭照护者及居家康复患者
Travel & Emergency Hygiene	Travellers, campers and people without bathing access	旅行及应急清洁	旅行、露营及不便洗浴人群

DIRECTIONS FOR USE

使用说明

01

Open the pack and remove one shampoo cap.
打开包装，取出一个浴帽。

02

Place the cap on the head and gently massage the hair and scalp until thoroughly cleansed.
将浴帽戴在头上，轻轻揉搓头发和头皮，直至清洁完成。

03

No rinsing is required; towel dry if necessary.
无需冲洗；如有需要，可用毛巾擦干。

04

Dispose of properly after use.
使用后请妥善丢弃。

01



02



03



04



“For a more comfortable care experience, the unopened pack may be warmed in a microwave.
为获得更舒适的护理体验，可将未开封产品放入微波炉中加热。”

”

01

Microwave power should not exceed 550W;
微波炉功率不得超过 550W；

02

2.Heating time should not exceed 15 seconds;
加热时间不得超过 15 秒；

03

Always check the temperature before use;
使用前请先测试温度，确认温度适宜；

04

Do not overheat or reheat, as this may cause burns.
请勿过度加热或重复加热，以免烫伤。

OTHER VINNER PRODUCTS

维尼健康其他产品

Body Care

身体护理系列

EN		CN	
Product	Code	产品	型号
Bathing Wipes	VLBBW8	产品名称	VLBBW8
Bathing Gloves	VLBBG8	免冲洗沐浴手套	VLBBG8
Rinse-Free Shampoo Cap	VLSHMC1	免洗浴帽	VLSHMC1
Continence Care Wipes	VLCONCW25	失禁护理湿巾	VLCONCW25
Pure Water Wipes	HKPWW60	纯水婴儿湿巾	HKPWW60

Body Decolonization

身体去殖系列

EN		CN	
Product	Code	产品	型号
CHG Wash Cloths	VLCWC	CHG 洗必泰身体擦拭巾	VLCWC
CHG Wash Gloves	Coming Soon	CHG 洗必泰身体擦拭手套	即将上市
CHG Shampoo Cap	Coming Soon	CHG 洗必泰免洗浴帽	即将上市

Clinical Disinfection

医疗消毒系列

EN		CN	
Product	Code	产品	型号
Surface Disinfectant Wipes	VLCDW	季铵盐器械消毒湿巾	VLCDW
Alcohol Disinfectant Wipes	VLADW	酒精器械消毒湿巾	VLADW
Peracetic Acid Wipes	Coming Soon	过氧乙酸消毒湿巾	即将上市
Hydrogen Peroxide Disinfectant Wipes	Coming Soon	过氧化氢湿巾	即将上市

OEM/ODM

代工生产/定制代工

EN			CN		
Product Category	Customizable Products	Formula Availability	产品类别	可定制产品	配方研发
Baby & Mother Care Wipes	Baby wipes, maternity wipes, breast cleansing wipes, sensitive skin wipes	✓	母婴护理湿巾	婴儿湿巾、产妇护理湿巾、乳房清洁湿巾、敏感肌湿巾	✓
Personal Hygiene Wipes	Feminine hygiene wipes, toilet wipes, deodorizing wipes, travel wipes	✓	个人卫生湿巾	女性护理湿巾、厕用湿巾、除味湿巾、旅行湿巾	✓
Medical & Healthcare Wipes	Patient cleansing wipes, pre-op cleansing wipes, wound surrounding skin cleansing wipes, medical care wipes	✓	医疗及健康护理湿巾	患者清洁湿巾、术前清洁湿巾、伤口周围皮肤清洁湿巾、医疗护理湿巾	✓
Disinfection Wipes	Surface disinfectant wipes, alcohol wipes, non-alcohol disinfectant wipes, sporicidal wipes	✓	消毒湿巾	表面消毒湿巾、酒精湿巾、无酒精消毒湿巾、孢子杀灭型消毒湿巾	✓
Pet Care Wipes	Pet cleansing wipes, deodorizing wipes, paw wipes, ear and eye wipes	✓	宠物护理湿巾	宠物清洁湿巾、宠物除味湿巾、宠物足部湿巾、耳眼护理湿巾	✓
Household & Cleaning Wipes	Kitchen wipes, bathroom wipes, floor wipes, glass wipes, multi-purpose cleaning wipes	✓	家居清洁湿巾	厨房湿巾、浴室湿巾、地板湿巾、玻璃湿巾、多用途清洁湿巾	✓
Industrial & Specialty Wipes	Equipment cleaning wipes, heavy-duty cleaning wipes, anti-static wipes, customized functional wipes	✓	工业及特殊用途湿巾	设备清洁湿巾、重油污清洁湿巾、防静电湿巾、定制功能湿巾	✓

VINNERLOVE provides a complete patient hygiene and infection prevention portfolio, covering daily body care, body decolonization, and clinical surface disinfection. From bedside bathing to professional healthcare disinfection, we support safer, cleaner, and more comfortable care environments.
VINNERLOVE 提供完整的患者清洁护理与感染预防产品组合，覆盖日常身体护理、身体去殖护理及医疗表面消毒。从床旁沐浴到专业医疗消毒，我们致力于帮助医疗及护理机构打造更安全、更洁净、更舒适的护理环境。

Vinner (Shenzhen) Health Products Limited.
维尼健康（深圳）股份有限公司
Room 608, Baoshan Times building, Shenzhen Longhua District, China. 518129
深圳市龙华区宝山时代大厦 608室

www.vinner-wipes.com
hongzihao@vinnercare.cn
liulizan@vinnercare.cn
TEL: +86 0755-61695825

